
# Minilight

Installation and use manual



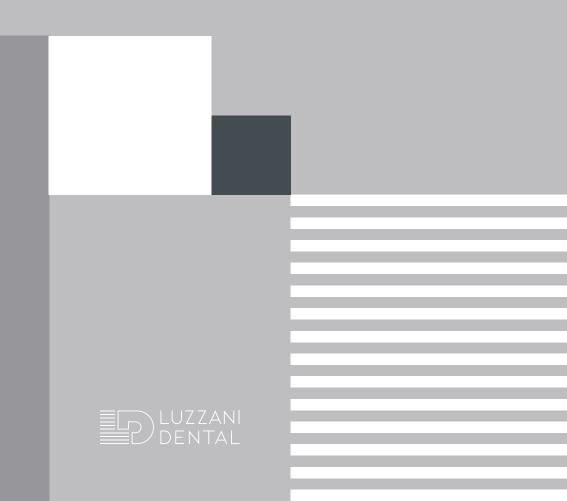

LUZZANI DENTAL S.R.L.

Via Torino, 3 20030 Senago (MI) Italy

www.luzzani.eu

E-mail: info@luzzani.com

Tel.: +39 029988433 Fax: +39 0299010379

ED. 5 REV. 0



### **CONTENTS**

0 — SYRINGE LEGEND
1 — WARNINGS
2 — PRODUCT DESCRIPTION
<ul> <li>2.1 — GENERAL</li> <li>2.2 — GENERAL CHARACTERISTICS</li> <li>2.3 — MODELS</li> <li>2.4 — CONTROLS</li> <li>2.5 — CE MARKING</li> </ul>
3 — IDENTIFICATION DATA AND WARRANTY
3.1 — MANUFACTURING BATCH 3.2 — WARRANTY
4 — PACKAGING
5 — TECHNICAL CHARACTERISTICS
5.1 — OPERATING CONDITIONS 5.2 — TRANSPORT AND STORAGE

### 8 — CLEANING AND STERILIZATION

9 — MAINTENANCE

10 — DISPOSAL AND SCRAPPING

11 — INFORMATION FOR THE DENTIST

12 — FAULT REPORT FORM

13 — SYMBOLS

14 — WIRING DIAGRAM

# 6 — INSTALLATION AND CONNECTIONS

ABSORBER CORD

**CONDITIONS** 

- 6.1 CONNECTION TO ELECTRICAL SYSTEM
  6.2 CONNECTION TO HYDRAULIC SYSTEM
  6.3 CONNECTION TO COMPRESSED

  AIR SYSTEM
  6.4 CONNECTION OF STRESS
- 6.5 NOTES FOR CORRECT CONNECTION

### 7 — NORMAL USE

7.1 — INSUFFLATION OF COLD WATER

7.2 — INSUFFLATION OF COLD AIR

7.3 — COMBINED INSUFFLATION
OF COLD WATER AND AIR (SPRAY)

7.4 — INSUFFLATION OF WARM WATER (present in version 6f and L)


7.5 — INSUFFLATION OF WARM AIR (present in version 6f and L)


7.6 — COMBINED INSUFFLATION OF WARM
WATER AND AIR (SPRAY)
(present in version 6f and L)

7.7 — LIGHT FUNCTION

7.8 — FIRST TIME USE AND USE AFTER LONG PERIODS OF NON-USE

### 0 —— SYRINGE LEGEND

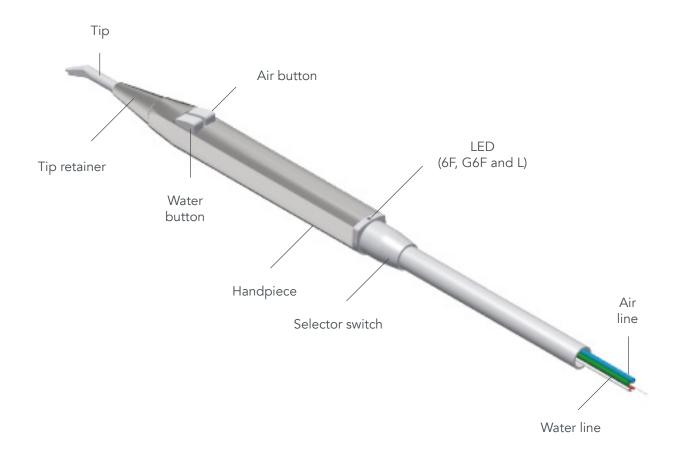

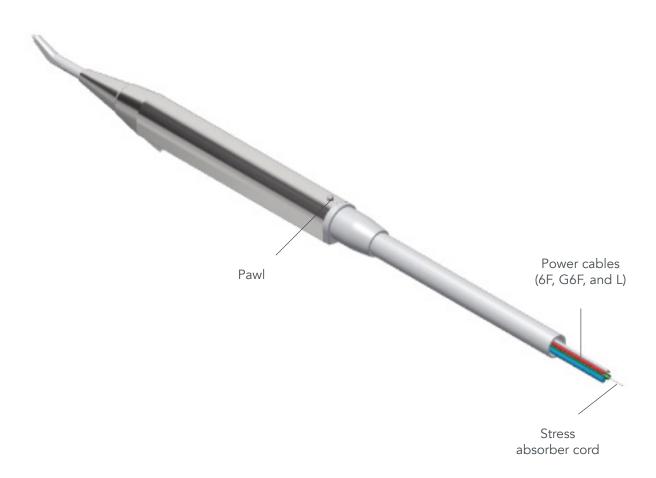

### 1 — WARNINGS

- 1.1 Any unauthorised tampering, modification or improper use immediately terminates the warranty and exonerates our company from all liability for injury or damage to persons, animals or property that may be caused by such interference.
- 1.2 To ensure maximum user and patient safety, the transformer used must be rated Safe Extra Low Voltage (SELV) with CE marking for medical use, in compliance with IEC 80.601 and IEC 60.601.
- 1.3 Connection to terminals with output voltage exceeding 24 V could irreparably damage the device.
- 1.4 The stress absorber cord must be anchored to the dental unit. This cord is designed to absorb any strains, thus preventing abnormal traction on the electrical or air/water line connections.
- 1.5 Before use, the electric, water and air lines must be correctly connected. The syringe could be damaged if used when not connected to the water and air supplies. The lines must be connected correctly since inverting them would damage the syringe.
- 1.6 To meet the requirements of Directive 93/42/EEC, the company has established a post-marketing surveillance procedure to monitor any problems generated by the use of our products. The attached form allows you to report any faults and suggest improvements which will be considered for subsequent versions of the product.
- 1.7 With each syringe, the package also includes a User's manual which also includes a "Fault Report Form". Since this is required by law, the user must receive this User's manual. Therefore, the syringe installer is responsible for delivering this User's manual to the dentist. Directive 93/42/EEC requires product traceability: therefore, our customers are asked to ensure that, in case of emergency, we can identify the end customer to whom the product has been sold.
- 1.8 Manufacturers and/or installers of dental units are required to comply with all the prescriptions outlined in this document.
- 1.9 Use the Minilight syringe only for the applications described in the instructions for use.
- 1.10 This product must only be installed by qualified persons.
- 1.11 Never modify the syringe in any way. This is strictly forbidden.
- 1.12 Use only original parts produced or approved by the manufacturer. If non-OEM accessories or consumables are used, the company cannot guarantee safe operation and function. No claims can be made for damages resulting from improper use.
- 1.13 Disposable tip adapters are not part of the Minilight Syringe and may compromise proper functioning.

- 1.14 Do not use the device in close contact with anaesthetic gases or in highly oxygenated environments (with an oxygen content >25%) or in areas where there is a risk of explosion.
- 1.15 Do not perform any maintenance procedures not indicated in the manual.
- 1.16 Before using the syringe, make certain that the water and air supplies have been activated.
- 1.17 National regulations regarding dental unit water and air quality must be met.
- 1.18 The air used must be dry, clean and free of oil.
- 1.19 The Minilight syringe meets the requirements laid out in European Council Directive 93/42/EEC: Class II A
- 1.20 The Minilight syringe must not be used near or set on other devices.
- 1.21 Use only accessories, cables, transducers specified or supplied by Luzzani Dental.
- 1.22 Do not use mobile RF communication devices within a distance less than 30 cm from the device.

### 2 — PRODUCT DESCRIPTION

### 2.1 — GENERAL

The Minilight syringe is a medical device designed to blow air and water (separately or together, at room temperature or at body temperature) to clean and/or dry the oral cavity during any dental procedure. It has been designed for use in dentist offices and dental clinics and is built into dental unit used exclusively by dentists. Product life — under proper maintenance conditions — is 5 years.

### 2.2 — GENERAL CHARACTERISTICS

- The Minilight syringe is a medical device for dental use (class II a).
- Protection rating against direct contacts: B
- Temporary operation: 10 sec. ON/ 20 sec. OFF.

The Minilight syringe has been designed using the latest ergonomic concepts for easy use and immediate cleaning and sterilisation. Both the tip and handpiece can easily be removed for perfect autoclave disinfection and sterilisation at 134°C (see point 8). Several handpieces of different shape are available. Choose the shape you need: angled or stylet. The air and water can also be heated to body temperature, thus sparing the patient even the slightest discomfort induced by insufflation of products at ambient temperature.

### 2.3 — MODELS

The models are differentiated into versions based on the

number of functions provided:

- 3F: cold water/air/spray
- 6F: cold water/air/spray warm water/air/spray
- G3F: cold water/air/spray (air and water inverted)
- G6F: warm water/air/spray (air and water inverted)
- A: cold water or air only
- L: cold/warm water, air and spray + light

#### Note:

Models 3F, G3F, 6F, G6F and L can be recognized by the printing on the back of the syringe.



Various handpieces, that differ in shape, can be mounted on all versions of the Minilight syringe.

The shape of the handpiece can be:

### **STYLET**



**ANGLED** 



This difference lets the dentist choose the tool ergonomically best suited for the purpose. To guarantee maximum hygiene and atoxicity, the handpieces are made of stainless steel. The devices are produced entirely in our workshop, with a tested, constantly updated work cycle using the most sophisticated machinery compliant with current quality system directives (UNI EN ISO 13485 certified).

### 2.4 — CONTROLS

Every syringe and all of its parts undergo duly documented, 100% complete functional and safety testing to ensure that the technical and functional design requirements are fully met.

### 2.5 — CE MARKING

All products bear CE marking both on syringe handpiece and inside. (batch number, autoclave symbol, Luzzani Dental logo, product name, CE marking with Notified Body number). The User's manual supplied with the product also includes details of our company, the main product characteristics and instructions for correct use and maintenance.

# 3 — IDENTIFICATION DATA AND WARRANTY

### 3.1 — MANUFACTURING BATCH

A number, marked on the inside of each product, identifies the production batch; the number is printed on the central body of the syringe. This number uniquely identifies the production batch thus always guaranteeing traceability of the product and each of its components, with relative test sheets.

### 3.2 — WARRANTY

The product is guaranteed by our company for 12 months from the date of the delivery document. The warranty covers any device manufacturing defects (materials) and is limited solely to the replacement of defective parts, performed in our workshop. The product must be sent to our premises at the expense of the customer. For the warranty to be valid, the product must be returned intact, complete and showing with no signs of tampering. The syringe has no functional expiration date; its expected life span is 5 years.

### 4 — PACKAGING

The product is shipped in suitable packaging to prevent problems during transport. The packaging consists of a plastic bag containing the Minilight syringe. A protective handpiece sheath and tip-saving tube protect the syringe during transport. Several bags are placed in one box. The Minilight syringe comes ready for connection to the dental unit once all packaging has been removed. For the correct use of the syringe, the protective tube must also be removed from the tip.

## À

### IMPORTANT NOTE:

With each syringe, the package also includes a User's manual which also includes a "Fault Report Form". Since this is required by law, the user must receive this User's manual. Therefore, the syringe installer is responsible for delivering these forms to the dentist. Directive 93/42/EEC requires product traceability: therefore, our customers are asked to ensure that we can identify the end customer to whom the product has been sold.

### 5 — TECHNICAL CHARACTERISTICS

FUNCTIONS		6F	3F	L
Supply voltage	V~	24	***	24
LED voltage	Vdc	***	***	3,5
Max. absorbed current	А	4,3	***	4,3
Electrical power	W	103	***	103
Water supply pressure	Кра	250	250	250
Air supply pressure	Кра	450	450	450
Water flow rate	Cc/mm	110	110	110
Air flow rate	NI/mm	10	10	10
International protection		IP40	IP40	IP40

### 5.1 — OPERATING CONDITIONS

Ambient temperature	10°C / +45°C
Relative air humidity	30% / 85%
Atmospheric pressure	80 Kpa - 106 Kpa

### 5.2 — TRANSPORT AND STORAGE CONDITIONS

Temperature	-20°C / +60°C
Relative Humidity	30% / 85%
Atmospheric Pressure	50 Кра - 106 Кра

# 6 — INSTALLATION AND CONNECTIONS

### 6.1 — CONNECTION TO ELECTRICAL SYSTEM

The connection involves hooking up the two electrical terminals to the transformer, arranged to provide an output of 24 V (see wiring diagram 13.1). Model L only: voltage for LED light 3.3 Vdc: black wire (+), white wire (-)



### NOTE

- The electric power supply and transformer used must be rated Safe Extra Low Voltage (SELV) for medical use, in compliance with the mandatory regulations.
- Connection to terminals with output voltage exceeding 24
   V could irreparably damage the unit and compromise safety.

### 6.2 — CONNECTION TO HYDRAULIC SYSTEM

The syringe's green line must be hooked up to the water supply.



### NOTE

- The operating pressure is 250 kPa. A slight increase in pressure would only increase the power of the jet, but would not create any problems and is not dangerous.
- The water used must be potable water, filtered (<25  $\mu$ m) and free of bacteria, etc.
- For the syringe to function properly, the water pressure must not be lower than indicated.

### 6.2 — CONNECTION TO COMPRESSED AIR SYSTEM

The syringe's blue line must be connected to the compressed air system.



### NOTE

- The recommended operating pressure is around 450 kPa.
- When using the Minilight syringe, national regulations regarding water and air quality must also be met.
- The air must be medical grade, dry and free of oil and bacteria a 5µm air filter is recommended.

### 6.4 — CONNECTION OF STRESS ABSORBER CORD

The stress absorber cord must be anchored to the dental

unit. This cord is designed to absorb any strains, thus preventing abnormal traction on the electrical or air/water line connections. The manufacturer cannot be held liable for malfunctions caused by failure to anchor the stress absorber cord.


### 6.5 — NOTES FOR CORRECT CONNECTION

- Before carrying out functional tests, both the electric, water and air lines must be correctly connected.
- Using the syringe without connecting the water and air supplies could damage the syringe.
- The lines must be connected carefully since inverting them could damage the syringe.

### 7 — NORMAL USE


### 7.1 — INSUFFLATION OF COLD WATER

To blow cold water into the operating field, just press the left button on the handpiece, symbol:

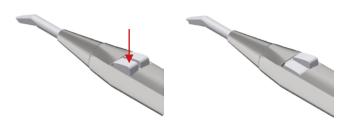


### 7.2 — INSUFFLATION OF COLD AIR

To insufflate cold air into the operating field, just press the right button on the handpiece, symbol:

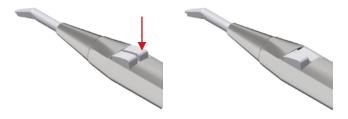

# 7.3 — COMBINED INSUFFLATION OF COLD WATER AND AIR (SPRAY)

To blow a combination of cold air and water (spray), press both buttons on the handpiece at the same time:

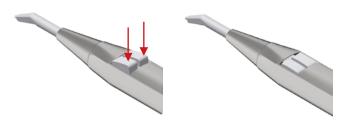


7.4 — INSUFFLATION OF WARM WATER (present in version 6f and L)

To blow warm water into the operating field, turn the selector switch at the base of the handpiece to the right (the green

### LED lights up)



and press the left button on the handpiece:

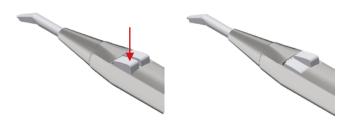

7.5 — INSUFFLATION OF WARM AIR (present in version 6f and L)

To blow warm air into the operating field, turn the selector switch at the base of the handpiece to the right (the green LED lights up):



and press the right button on the handpiece:




7.6 — COMBINED INSUFFLATION OF WARM WATER AND AIR (SPRAY) (present in version 6f and L)

To blow a combination of warm air and water (spray), turn the selector switch at the base of the handpiece to the right (the green LED lights up)



and press both buttons on the handpiece at the same time:

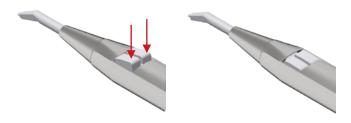

### 7.7 — LIGHT FUNCTION (present in version L)

This function is present in the Minilight L version. In this model, the handpiece has an optical fibre that conducts the light generated by a LED located in the body of the syringe. We recommend using a switch to turn the LED on and off. For connections, see wiring diagram 13.2.



#### NOTE

- The sole function of the selector switch is to preselect warm operations. The water and/or air are only heated at the moment in which they are actually used.
- The cleaning (or line washing) procedures must always be performed in the cold position.

### WARNING

Do not use the tip improperly. Remove and sterilise the tip after each patient.

### **IMPORTANT**

Air and water must be able to flow freely from the tip. Do not rest the tip on the tooth or on an object. Do not press the tip against impression materials as they could cause obstruction

### 7.8 — FIRST TIME USE AND USE AFTER LONG INTERVALS



- Sterilize the handpiece and all accessories before
  use
- After prolonged periods of inactivity, clean, treat and sterilise the handpiece.

# BEFORE EACH PATIENT

- 1. Make certain the handpiece has been sterilised.
- 2. Adjust the supply of fluids from the dental unit (see table in point 5).
- 3. Press the air button and make certain that there is a clearly perceptible jet of air.
- 4. Check the water flow rate.
- 5. Use only filtered water that is free of oil and microorganisms.

6. Check the tip for any obstructions or deposits. Clean if necessary.

#### NOTE

- Flush out the syringe at the beginning of each work day (minimum flushing time: 2 minutes) and before each patient (minimum flushing time: 20-30 sec.).
- Immediately upstream of the syringe, install filters able to retain the microorganisms coming from the hydropneumatic circuit.

### 8 — CLEANING AND STERILIZATION



After each use on a patient, the handpiece and tip of the syringe MUST be cleaned and sterilised to guarantee maximum hygiene.

Sterilisation symbol on the handpiece:



To do this, proceed as follows:

Disconnect the tip by unscrewing the tip retainer


1. unscrew the tip retainer



2. withdraw the tip retainer




3. withdraw the tip

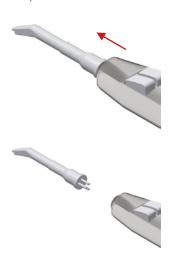

or remove the entire handpiece by pressing the button on the lower part of the handpiece and pulling upwards.

1. press the pawl to release the handpiece

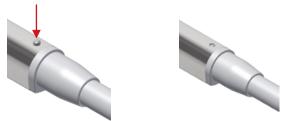

### 2. withdraw handpiece



wipe with a damp cloth, removing any stains. Set in a steam autoclave at 134°C for AT LEAST 3 minutes (in compliance with CEI EN 13060).

### A — WARNINGS



The syringe should always be sterilised, even before using it for the first time.

- Inappropriate sterilisation is hazardous for the patient and for the operators.
- Do not perform spray disinfection.
- Do not immerse in disinfectant liquids.
- Do not perform cold or hot air chemical sterilisation.
- The personnel performing the task must be skilled and specially trained.
- Use disinfectant according to the specifications on the manufacturer label.
- Do not use chlorine-based liquids.
- When simultaneously sterilising more than one item in an autoclave, check that the load does not exceed the maximum allowed.

### B — PREPARATION

Eliminate surface dirt using a disposable paper towel. Clean the inside of the lines by running air and water through the syringe for about 30 seconds. Remove the stainless-steel handpiece by pressing the button on its terminal section. Unscrew the tip retainer and remove the tip.

### C — MANUAL CLEANING

Use a disposable paper towel and potable water to remove any impurities or dirt that may be present

### D — AUTOMATIC CLEANING 「本」



Not envisaged

E — MANUAL DISINFECTION

product technical data sheet).

Recommended disinfectants:

- FD 322 Durr
- Mikrozid AF liquid

### F — AUTOMATIC DISINFECTION



Not envisaged

### G — MANUAL DRYING

Dry with disposable paper towelling. Dry with clean, dry, uncontaminated compressed air, inside and out, continuing until completely dry. Do not dry with hot air.

Disinfect only with a disposable cloth and the permitted disinfectant (following the instructions on the label and

### H — AUTOMATIC DRYING

Not envisaged

### I — MAINTENANCE AND CONTROL

No special maintenance is necessary. There is no objective period of time that limits the useful life of the handpiece: visually check for damage and signs of wear, and if found, replace the part.

### L — PACKAGING

Use heat sealable film-paper sterilisation pouches of appropriate size.





The handpiece and tip can be sterilised.

Sterilise in a class B steam autoclave in compliance with EN 13060 ISO 17665-1.

3-phase sterilisation with fractional vacuum system at 134°C +/- 1°C at a pressure of 2.13 bar, applying a 4-minute delay. Never exceed 134°C. The autoclave must be validated.

### N — STORAGE

No particular requirements apart from storage in the sealed, sterilised pouches. Store in a suitable place that is dry, out of direct sunlight and possibly with low bioburden.

### 9 — MAINTENANCE

The instrument requires no specific maintenance apart from normal cleaning and sterilisation as described in the previous paragraph.

### 10 — DISPOSAL AND SCRAPPING



The product does not contain dangerous or toxichazardous components. Separate waste collection is required for electrical equipment. Follow the regulations in force in your country.

### 11 — INFORMATION FOR THE DENTIST



The dental unit manufacturer is required to deliver the Luzzani Dental syringe User's manual to the end user.

### 12 — FAULT REPORT FORM

To meet the requirements of Directive 93/42/CEE as amended, the company has established a post-marketing surveillance procedure to monitor any problems generated by the use of our products. This commitment includes the requirement that both user and manufacturer inform the competent authorities of any incident caused to patient or user by malfunction or deterioration in the characteristics and/or performance of a device, as well as any inadequacy in the instructions for use. We kindly ask you to inform us of any anomalies by sending us the sheet attached to the last page of this manual.

### 13 — SYMBOLS

### **SYMBOLS**



Do not overturn



Fragile



Keep dry



Type B device



Alternating current



General warnings



Separate collection for electrical and electronic equipment



Manufacturer



Double insulation



Batch number




Sterilise




Consult User's manual

### 14 — WIRING DIAGRAM

### 14.1 – VERSION 6F

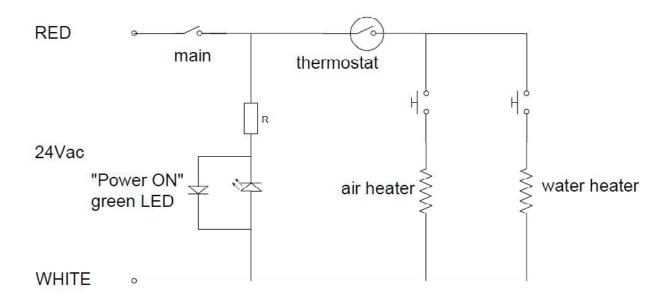

### 14.2 – VERSION L

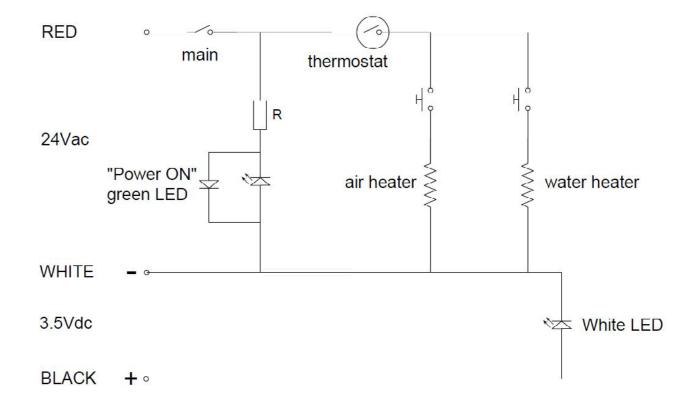

### FAULT REPORT FORM

PRODUCT	
TYPE B <sub>A</sub>	ATCH
REPORTED BY	
COMPANY	
TYPE OF REPORT	
ANOMALY SUGGESTIONS	
DESCRIPTION	
NOTES	
DATE	SIGNATURE

SEND TO:

LUZZANI DENTAL SRL Via Torino 3 - Senago (MI) - ITALY Tel. +39 02 99010379)