



SHIELD FORCE PLUS



Fast pain relief


The reaction of the 3D-SR-Monomer provides relief in seconds.
Block 1 - Closure of the tubules.

Permanent pain relief

The resulting protective layer on the tooth surface is light-cured.
Block 2 - Sealing and permanent protection for up to 3 years.

Handling

Extremely easy to use. No rinsing or massaging. Simply apply, leave on for 10 seconds, air blow and light cure. Done!

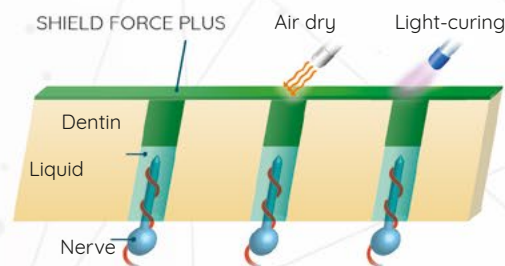


Find out more:



Thermal cycle test (10,000 cycles) at 4°C-60°C **

The sealing layer of SHIELD FORCE PLUS remains intact after the thermal cycle test.

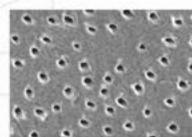

Hypersensitive dentin model




After application of SHIELD FORCE PLUS




After the thermal cycle test (10,000 x)


Sealing layer of the Gluma* Desensitizer (Kulzer) washes out during the thermal cycle test.

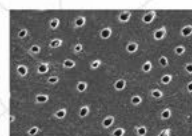

Hypersensitive dentin model

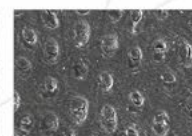

After the application of Gluma* Desensitizer

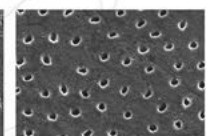

After thermal cycle test

* Not an own trademark of TOKUYAMA Dental Corp.
** Source: TOKUYAMA DENTAL R&D



Keyfacts Desensitiser

- > Quick and easy handling
- > Excellent sealing effect
- > Durable & resistant protective coating