

TF4 | RECIPROCATING SYSTEM

Material : Niti Gold standard Length : 21mm,25mm,31mm

Rotary Speed : 350 rpm Torque : 2.0~3.0 N/cm

Sterilization: ≤126°C

Package : 4pcs blister/box assorted/single size

● 20/.07 ● 25/.07 ● 35/.06 ● 45/.05

Cross Section:



Guidelines For Use

- If a 10 K-file is very resistant to movement, use #20.
- If a 10 K-file moves to length easily, is loose or very loose, use #25.
- If a 20 hand file or larger goes to length, use #40





TF3 and TF4

Material : Niti Gold standard Length : 21mm,25mm,31mm
 Rotary Speed : 350 rpm Torque : 2.0~3.0 N/cm
 Package : 3 or 6 pcs per box/blister assorted/single size

● 20/.06 ● 25/.08 ● 40/.08

Cross Section:




What is the difference between TF3 and TF4?




The single-file simplicity of TF3 is enhanced in TF4. Made with advanced metallurgy.

TF4 File is 50% more resistant to cyclic fatigue than TF3

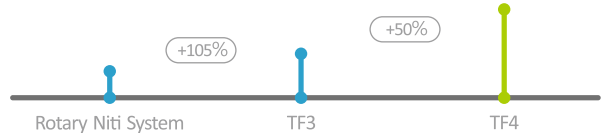

Enhanced file flexibility thanks to the Gold technology

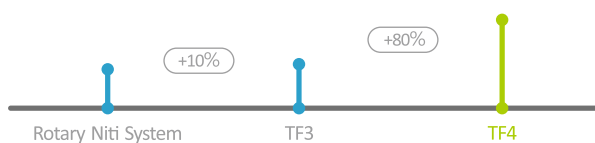

Shaping working time(Takes the cutting efficiency to a higher level)

