



DuoTEMP[®] 
dual cure temporary filling material

Double security keeps you smiling!

- Dual-cure
- High marginal integrity
- Easy removal
- One material with multiple indications

 **COLTENE**

DuoTEMP® is a dual-cured, eugenol-free, white temporary filling material formulated from zinc oxide/zinc sulphate for optimal protection from the infiltration of bacteria. DuoTEMP is one of the first all-in-one temporary filling materials that can be cured using both light and saliva. This solution delivers greater flexibility, allowing the patient to immediately resume normal daily activities. DuoTEMP also contains a pleasant flavor that keeps the dentist and patient smiling.

► **Double security**

- **Dual-cure**
After only 40 s of light curing, the patient can immediately load the material. Afterwards, DuoTemp will continue to completely self-cure on its own with saliva.
- High marginal seal due to slight expansion
Reduces infiltration of bacteria

► **One material with multiple indications**

- Reduces inventory and minimizes confusion
- Temporary fillings of class I and II cavities
 - Temporary fillings for endodontic applications
 - Temporary fillings for inlay / onlay preparations

► **Natural tooth-like feeling**

DuoTEMP's high surface hardness feels like a natural tooth to minimize irritation and maximize patient comfort.

► **Additional benefits**

- Contains zinc oxide and fluoride for additional protection
- Eugenol-free – no effect on the curing of composites
- Radiopaque


Order Information:

5831 DuoTEMP Single Pack
1 syringe of 5 g

7090 DuoTEMP Eco Pack
5 syringes of 5 g

7235 Intraoral Separator, 5ml

► **Marginal Seal Comparison**

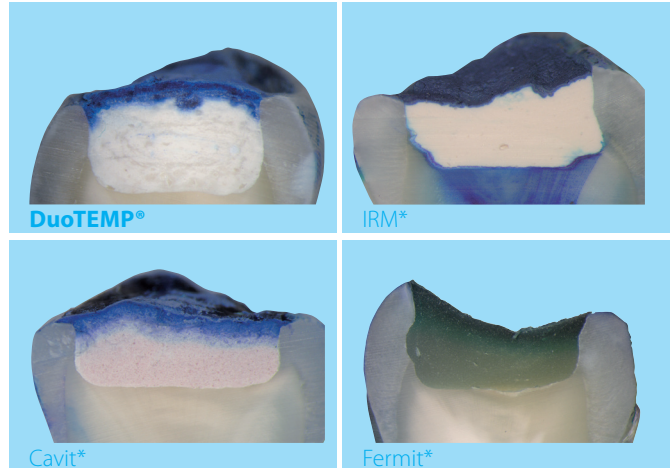

Cross-sectioned extracted human molars were filled with cured DuoTEMP and then immersed in methylene blue dye to determine the marginal sealing properties. DuoTEMP demonstrated excellent sealing properties.
Source: internal data

* Not a trademark of Coltène/Whaledent

► **Easy, time saving application**



One component without mixing. Good retention in the cavity. Easy removal.

Coltène/Whaledent Inc.
235 Ascot Parkway
Cuyahoga Falls, OH 44223 / USA
Tel. USA & Canada + 1 800 221 3046
+ 1 330 916 8800
Fax +1 330 645 8704
info.us@coltene.com

Coltène/Whaledent AG 
Feldwiesenstrasse 20
9450 Altstätten / Switzerland
Tel. +41 (0)71 757 53 00
Fax +41 (0)71 757 53 01
info.ch@coltene.com

Coltène/Whaledent GmbH + Co. KG
Raiffeisenstrasse 30
89129 Langenau / Germany
Tel. +49 (0)7345 805 0
Fax +49 (0)7345 805 201
info.de@coltene.com

coltene.com
COLTENE