

Temp-Bond[™]

Temporary dental cement – Automix syringes



Temporary attachment. *Permanent reliability*

- High bond strength
- Easy to remove
- Flexible delivery choices
- Flows and mixes easily



Together, we're more.[™]

Temp-Bond™


A history of innovation and success



With more than 50 years in the market, Temp-Bond is a brand that has its roots firmly planted in dental research. The original **Temp-Bond** launched in 1961 and was based on zinc oxide-eugenol chemistry, a classic setting reaction for temporary cementation with over 100 years of clinical history.




In 1989, Kerr introduced **Temp-Bond NE** as a solution for patients with eugenol allergy. The eugenol-free zinc oxide cement replaces the eugenol with carboxylic acids.

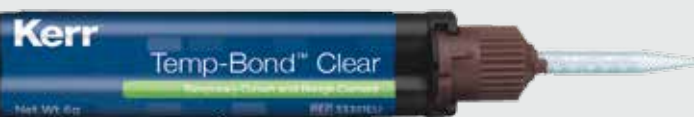

The latest instalment in these advances is **Temp-Bond Clear**. The formula introduces a unique redox initiator system and low refractive index fillers that give the cement a highly stable colour for a natural aesthetic look.

Material set times

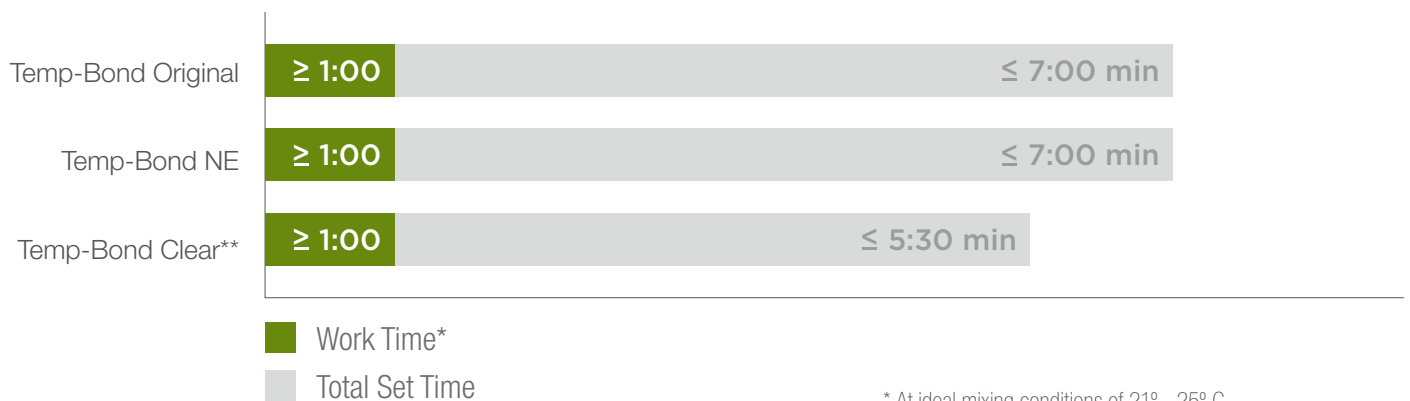

* At ideal mixing conditions of 21° - 25° C.
Work time listed for automix syringes only

** Temp-Bond Clear can be light cured for a more immediate set.

Features

High bond strength. Designed to prevent leakage.

Easy to remove. Minimises potential damage to the prep or temporary.

Flexible delivery choices. Temp-Bond Original, Temp-Bond Clear and Temp-Bond NE are available in automix syringe. Moreover, Temp-Bond and Temp-Bond NE are also available in Unidose® and tubes, offering flexible delivery options.

Flows and mixes easily. Delivers optimal consistency for solid, complete seating of restorations.

Proven performance you can trust

Throughout the years, the Temp-Bond family has affirmed itself as the leader in temporary cementation with its three products:

Temp-Bond – the temporary cement for trial cementing restorations or cementing temporary crowns and bridges.

Temp-Bond NE – maintains Temp-Bond's key characteristics while being non-eugenol therefore suitable for patients with eugenol allergy. It does not inhibit the polymerisation of permanent resin cements and acrylic temporaries.

Temp-Bond Clear – also eugenol free, it is the most translucent resin based provisional cement in the market for superior aesthetic results.

All your cementation needs taken care of

Effective cementation requires the right product depending on the indication, procedure or substrate you are addressing. Discover Kerr's leading permanent and temporary cements. Different solutions, same quality.

		Temp-Bond™ Temporary Cement	Maxcem Elite™ Self-Adhesive Resin Cement Permanent Cement	NX3 Universal Adhesive Resin Cement Permanent Cement
RECOMENDED INDICATIONS				
Permanent	CAD/CAM	-	Best	Best
	Inlays, Onlays and Crowns	-	Acceptable	Acceptable
	Veneers	-	Acceptable	Best
	Implant Crowns	-	Acceptable	Acceptable
	Posts	-	Best	Acceptable
	Maryland Bridges	-	Acceptable	Best
	Provisional Restoration	Best	-	-
	Temporary Crown and Bridges	Best	-	-
RECOMENDED SUBSTRATES				
Metal/PFMs			Best	Acceptable
Ceramic/Hybrid (VITA ENAMIC®)			Acceptable	Best
Zirconia (Bruxir®, Lava™ Zirconia)			Best	Acceptable
Lithium Disilicate Ceramics (IPS e.max®)LI			Best	Best
Composite (Lava™ Ultimate)			Acceptable	Best
All-ceramic (IPS empress)			Acceptable	Best

For a full list of indications and contraindications, please refer to the respective product Instructions for Use.

Temp-Bond™

Temporary Cement



Temp-Bond – Automix Syringe

CONTENTS: 2x 11.8g Syringe; 20x Mixing Tips.

Part No.	Description
33684	Temp-Bond – Automix Syringe

Temp-Bond – Refills

Part No.	Description
60254	Temp-Bond – 50g Tube Base
60255	Temp-Bond – 15g Tube Accelerator
61081	Temp-Bond – 13g Modifier

Temp-Bond – Introductory Pack

CONTENTS: 1x 50g Tube Base; 1x 15g Tube Accelerator; 1x 13g Tube Modifier; 1x Mixing Pad.

Part No.	Description
60058	Temp-Bond – Introductory Pack (with Modifier)

Temp-Bond – Standard Pack

CONTENTS: 1x 50g Tube Base; 1x 15g Tube Accelerator; 1x Mixing Pad.

Part No.	Description
61087	Temp-Bond – Standard Pack (without Modifier)

Temp-Bond – Unidose

CONTENTS: 50x 2.4g Unidose.

Part No.	Description
31377	Temp-Bond Unidose – 120g



Temp-Bond NE – Automix Syringe

CONTENTS: 2x 11.7g Syringe; 20x Mixing Tips.

Part No.	Description
33685	Temp-Bond NE – Automix Syringe

Temp-Bond™ NE

Non-Eugenol Temporary Cement

Temp-Bond NE – Unidose

CONTENTS: 50x 2.4g Unidose.

Part No.	Description
31973	Temp-Bond NE Unidose – 120g

Temp-Bond NE

CONTENTS: 1x 50 ml Tube Base; 1x 15ml Tube Accelerator; 1x Mixing Pad.

Part No.	Description
60256	Temp-Bond NE



Temp-Bond Clear – Accessories

Part No.	Description
33361	50x Automix Tip Refill

Temp-Bond™ Clear

Temporary Crown and Bridge Cement

Temp-Bond Clear – Automix Syringe

CONTENTS: 1x 6g Syringe; 10x Mixing Tips.

Part No.	Description
33351EU	Temp-Bond Clear – Automix Syringe